

# Staining Table FT 850-AK-ULS - Standing Workplace



This specialized work table provides a pollutant-free laboratory environment while staining histological slide preparations. Airflow travels over the entire working area with an integrated fan as well as a permanganate carbon filter and a HEPA filter. The air is vacuumed out downwards which provides optimal removal of unsanitary vapours. MAC (maximum allowable concentration) values are below target throughout the entire vacuum-effected area. Stable stainless steel construction with perforated working plate (alternatively available with openings for glass cuvettes at an extra charge). Sub construction with height-adjustable feet. 3-side plexiglas cover with stainless steel backside. Soundproofed silent exhaust fan with operation hour counter. Removable water container approx. 0.8 Liter.

## **Technical Data:**

Dimensions W/D/H :	850 x 850 x 1750 mm
Air speed:	min. 0,4 m/sec on the perforated working plate
Material:	Stainless steel (material no. 1.4301 – ANSI 304)
Waste drain:	DN 50 (optional)
Voltage:	230 Volt / 1 Ph / 50
Type:	FT-850-AK-ULS
Article No.:	060.260.150

*More information see back page*

**KUGEL medical GmbH & Co. KG**  
Hermann-Köhl-Str. 2A  
DE - 93049 Regensburg/Germany

Phone +49/941/20 86 48-0  
Fax +49/941/20 86 48-29  
[www.KUGEL-medical.de](http://www.KUGEL-medical.de)




# Staining Table FT 850-AK-ULS - Standing Workplace



Perforated working area with openings for glass cuvettes (customized solution)

## Accessories:

Water tap with long hose  
Type: EHMB-300  
Article no. 010.010.040

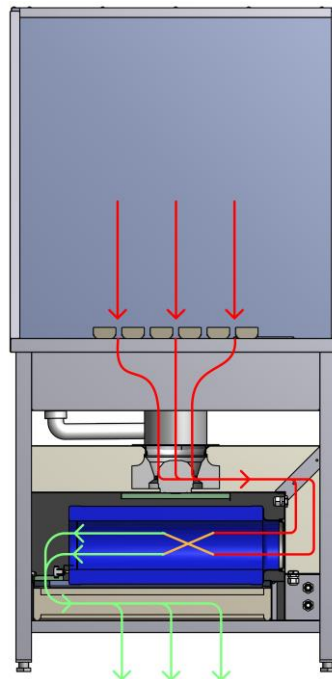

## Carbon filter (spare part)

Type: AKF-500-500/200  
Article No. 060.260.026

## HEPA Filter (spare part)

Type: HEPA H13-305-610  
Article no. 060.265.020

**KUGEL medical GmbH & Co. KG**  
Hermann-Köhl-Str. 2A  
DE - 93049 Regensburg/Germany

Phone +49/941/20 86 48-0  
Fax +49/941/20 86 48-29  
[www.KUGEL-medical.de](http://www.KUGEL-medical.de)

